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Before using the unit, please read the following instruction manual carefully. An improper use of the device and not as written in the manual can be dangerous for safety.

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The instrument is classified as electrical/electronic apparatus and must be disposed of accordingly. The unit is subject to waste separation and cannot be disposed of as urban waste under EEC directive 2002/96/CE. For more information, please contact the relative division of your local town council.

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**This unit must be used for laboratory applications only and for internal use only.**

The manufacturer declines all responsibility for any use of the unit that does not comply with these instructions.

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**This unit has been designed and manufactured in compliance with the following standards:**

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VELP reserves the right to modify the characteristics of its products to constantly improve their quality.

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The exclusive wireless data transmission from the RESPIROMETRIC Sensor to the Wireless DataBox™ can be set at selected intervals ranging from 30 minutes to 48 hours. The sampling time value, depends to the total duration of analysis, till 730 days. RESPIROMETRIC Sensor automatically stores the most recent values. Data transmission is ensured even when the Wireless DataBox™ is not connected to the PC. The Wireless DataBox™ can store unlimited results for up to 400 samples.

## 2 Assembling and installing

Upon receipt and after having removed the packaging, please check the integrity of the instrument. The box includes:

- Wireless DataBox™
- Power USB cable
- RESPIROSoft™
- Power supply USB

The Wireless DataBox™ comes already prepared and doesn't need any action for the assembling it. Using the USB cable is possible connect it to the PC or to the power supplier when analyses are in progress. In this way is possible record the dates inside the unit without use the computer.

On the USB-pen drive is possible to find the driver if required by Windows after the first connection to the PC.

## 3 RESPIROSoft™

### 3.1 Installation

General rules for installing the RESPIROSoft™ Software:

In order to install the RESPIROSoft™ Software correctly, administrator rights for the PC must be available. The PC on which the RESPIROSoft™ Software will be installed must not be overloaded with other software. This is to prevent unstable operation of the RESPIROSoft™ Software.

To install the RESPIROSoft™ Software:

1. Place the USB-pen Drive in the USB port.
2. Start the installation running the function **setup.exe** from the USB-pen drive.
3. The installation guide is displayed. Click **Next** and follow the instructions. The RESPIROSoft™ icon will be created on the screen automatically.
4. Connect the Wireless DataBox™ to the PC.
5. After a few seconds the PC acknowledges the presence of a new hardware and will automatically start the installation procedure of the drivers. When requested, specify the path of the **Driver** folder on the USB-pen drive and proceed with installation following the instructions that appear on the screen.
6. Re-boot the computer to put the settings provided by the installation procedure into effect.

### 3.2 Turning on

Before turning on the instrument, all the steps described in the previous paragraphs must be carried out:

1. Turn on the PC and the screen.
2. Check that the RESPIROSoft™ has been installed.
3. Connect the Wireless DataBox™ to the PC using the USB cable.
4. Launch the RESPIROSoft™ (shortest way > double-click on the appropriate icon on the screen) (*Fig. 1*).
5. The "User" level password is "velp". Enter the password and press OK (*Fig. 2*).
6. The main window will appear (*Fig. 3*).

Fig. 1

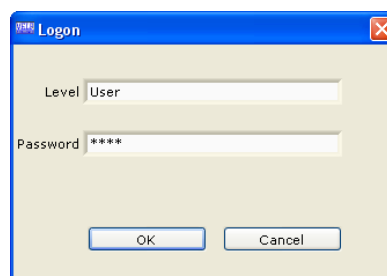
Fig. 2

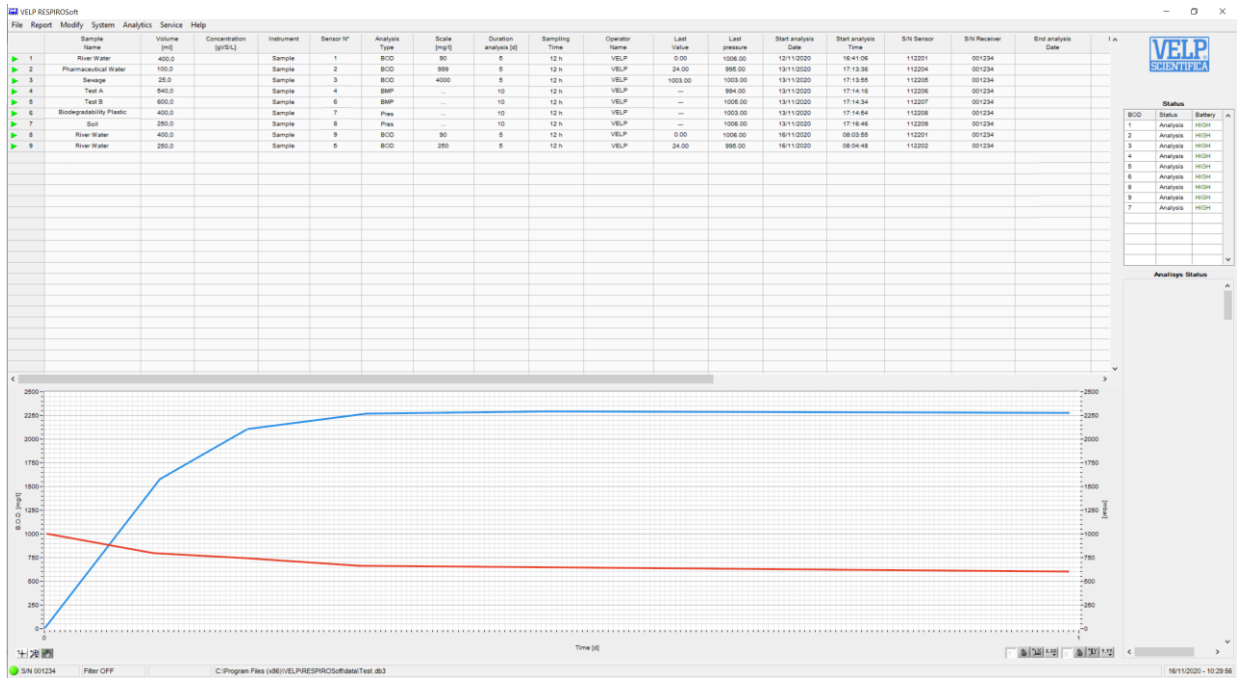
Fig. 3

## 4 RESPIROSoft™ Rapid Guide


This page shows an abstract of the RESPIROSoft™.

The purpose is to help the user identify the various possibilities offered by the instrument and how to reach them easily.

<b>Menu</b>	<b>Submenu</b>	<b>Description</b>
File	New	Create a new analysis data base and "Laboratory information" file.
	Open	Load an already existing database.
	Recent Files	Shortcut to recently opened databases.
	Save as	Save an analysis data base (including "Laboratory information" file) assigning a file name. If the file name already exists, it can be overwritten.
	Laboratory information	Enter the details of the laboratory concerned. This information will be included in the test report.
	Exit	Exit from the programme and return to the desktop.
Report	Create Report	Create a single or multiple test report with calculation of average and RSD values.
	Export Data	Export data of the rows selected to an .xls or .txt or csv file.
Modify	Filter	Create database filters.
	Standard column	Reset the order of the database columns to the default setting.
	Configure	Configure the sequence of the database columns, show/hide columns.
	Configure Copy	Allows to define fields which are copied with the Copy command
	Remove	Remove one or more rows from the database.
	Select all	Select all the rows of the database containing a value.
System	Languages	Change the language of the program.
Analytics	Standards	Display, modify, create or delete standards.
	Methods	Display, modify, create or delete analytical methods.
	BMP Calculator	Calculator for BMP Analysis
	Sampling configuration	Allows to create personalized sampling time and some other information such as analysis type and analysis duration
	Sensors DB	Allows to control name, serial number and type of already registered sensors (giving the possibility to remove them)
Service	Connect Databox	Manual search for Wireless DataBox™.
	Disconnect Databox	Disconnect the Wireless DataBox™ from the PC.
	Ermes	Enable/Disable ERMES connection
	Start AP	Start the AP mode on Wireless DataBox
	Acoustic signals	Select the operations that activate an acoustic signal.
	Login	Modify the software password and entry level.
	Diagnostic	Perform a diagnosis of the status of the instrument.
	Control panel	For VELP or authorized technical staff only.
Help	Software activation	Under construction.
	About	Contact VELP Scientifica directly by e-mail.
	Manual	PDF operating manual.

#### 4.1 Configure Copy

The menu Modify has the “Configure copy” option: the table allows to define fields which are copied with the Copy command. Fields which are copied are on the right, while fields which are not copied are on the left. Fields can be moved using the arrows buttons in the center.
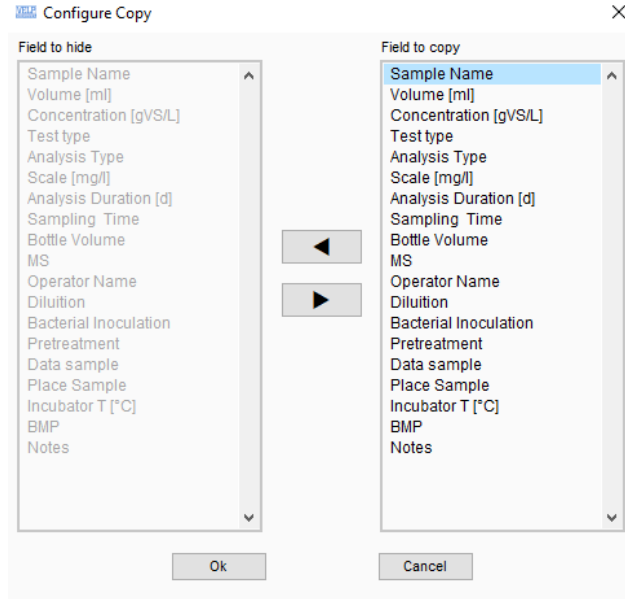
#### 4.2 Sampling Configuration

The menu Analysis has the “Sampling configuration” option: this table allows to add new custom sampling time and shows default sampling time (which can’t be modified).

In order to add a sampling time it’s necessary to login as super-user and fill the necessary columns: “test type” (BOD,BMP,Press), “analysis duration” in days and “sampling time” (in h or min). The “enable” column allows to show or hide the row with the previous information. The “name” column is optional.

Name	Enable	Test type	Analysis Duration [d]	Sampling Time	Time
	YES	BOD	5	6	h
	YES	BOD	5	12	h
	YES	BOD	7	6	h
	YES	BOD	7	12	h
	YES	BOD	14	6	h
	YES	BOD	14	12	h
	YES	BOD	21	12	h
	YES	BOD	21	24	h
	YES	BOD	28	12	h
	YES	BOD	28	24	h

#### 4.3 Sensors DB

The menu Analysis has the “Sensor DB” option: this table shows the association between the RESPIROMETRIC Sensor registration number and its S/N. It’s possible to remove an association right clicking and selecting the “Delete” option.

**Sensors Database**

Sensor Name	Sensor Serial	Sensor Type
1	16754935	BOD-BMP
2	1745891	BOD-BMP
3	4071645	BOD-BMP
4	6047	BOD-BMP
5	1746869	BOD-BMP
6	92475	BOD-BMP
7	560916	BOD-BMP
8	560918	BOD-BMP

## 5 DataBox™ Connection

In order to connect the Wireless DataBox™ to the PC, is necessary to use the USB cable, supplied with the unit.

1. After having connected it to the PC, select **Connect Databox** from the **Service** menu (Fig. 4).
2. Using the “Connect Receiver” menu with the Automatic Search option, it is possible to find and connect automatically the Wireless DataBox™ (Fig. 5).
3. When the unit is connected to the RESPIROSoft™, the software shows the S/N and a green led on the left bottom (Fig. 6).

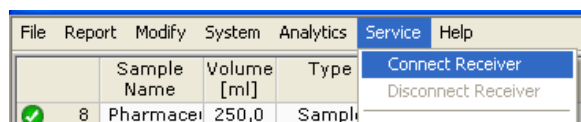
Fig. 4

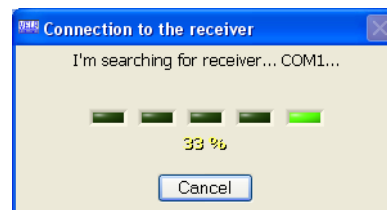
Fig. 5

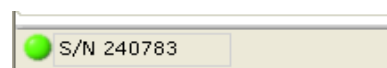
Fig. 6

## 6 RESPIROSoft™ Workflow

### 6.1 Preparing for the analysis

About the sample preparation, please see the analytical methods section of RESPIROMETRIC System Operating Manual. In the Analytics/Methods menu, it is possible to see a summary of methods reported in the analytical methods section of RESPIROMETRIC System Operating Manual.

#### 6.1.1 Filling the database

Several databases can be created and stored in the PC.

Each database is able to store a great number of analyses (if the speed is too low it is suggested to create a new database).

Select **File/New** to create a new database or select **File/open** to load an already existing database.

Each database is connected to a “Laboratory information” file (Fig. 7).

The information in this file will be included in the test reports.

The company logo can also be uploaded in order to customize the test reports.

The “Laboratory Information” file can be modified as many times as necessary.

Before starting the analysis, some of the columns must be filled in.

The following database information is mandatory and **must be entered** before starting the analysis:

<b>Sample Name</b>	Name of sample (alphanumeric field)
<b>Test Type</b>	Type of analysis to be performed
<b>Sensor N°</b>	Number used for identify the Sensor
<b>Analysis Type</b>	BOD, BMP or PRESSURE
<b>Scale</b>	Scale used for the analysis (just for BOD)
<b>Analysis Duration</b>	Total days of analysis
<b>Sampling Time</b>	Sampling time for the analysis

Each information will be explained in detail in the following paragraphs.

Once the analysis has been carried out the content of the database columns cannot be modified. All the other database columns are filled automatically by the instrument at the end of analysis. The information can be entered once and copied using **Copy/Add** (mouse right button) to add new row. The new row will be added in sequence. The submenu Configure Copy can be used to configure the parameters which are copied with **Copy/Add**.

In case of BMP, also Volume, Vbottle, MS and Incubator T must be compiled if the BMP value is required.

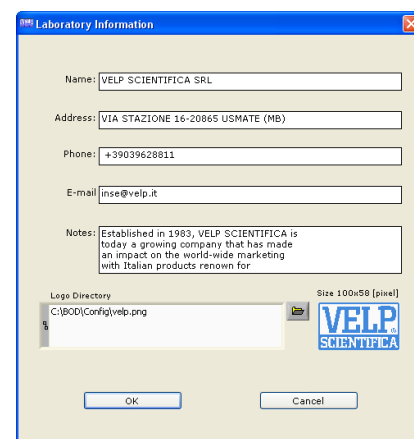
Fig. 7

### 6.1.2 Choosing the sample “Test Type”

The “Test Type” column offers the following options (Fig. 8):

<b>Sample</b>	Analysis on a sample.
<b>Standard</b>	When is necessary to compare the last obtained value with a range of results, e.g. with a standard in order to verify the correct working of the system

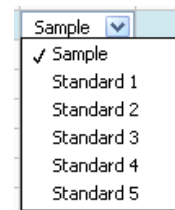
Fig. 8

In **Analytics/Standards**, it is possible to find some standard that are suggested for the test.

### 6.1.3 Choosing the “Sensor N°”

To use a Sensor with RESPIROSoft™, is necessary to register it with a number that is used to identify the unit. In order to register it, it is necessary to keep pressed “SET” on the RESPIROMETRIC Sensor until the word SET appears on the sensor. The display will then show the following screen (Fig. 9).

Give a number to the Sensor and press OK.

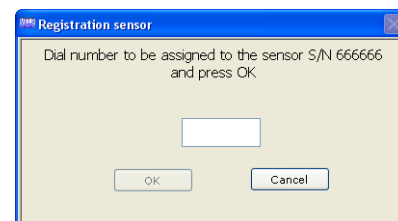
Fig. 9

In the column **Sensor N°** it is possible to select the units that have been already registered.

The registration of sensors is necessary only when use for the first time a RESPIROMETRIC Sensor with RESPIROSoft™. Otherwise, if the unit was already used, the number is already registered.

### 6.1.4 Choosing the “Analysis Type”

The “Analysis Type” column shows the following options:

<b>BOD</b>	Biochemical Oxygen Demand determination	BOD and Press graphs
<b>BMP</b>	Biochemical Methane Potential determination	BMP graph
<b>Press</b>	Pressure determination	Press graph

### 6.1.5 Choosing the “Scale”

Choose the **Scale** in case of BOD analysis, with the following options (Fig. 10):

Scale	Volume of sample
A: 0 ÷ 4000 mg O <sub>2</sub> /l	25 ml
B: 0 ÷ 999 mg O <sub>2</sub> /l	100 ml
C: 0 ÷ 600 mg O <sub>2</sub> /l	150 ml
D: 0 ÷ 250 mg O <sub>2</sub> /l	250 ml
E: 0 ÷ 90 mg O <sub>2</sub> /l	400 ml

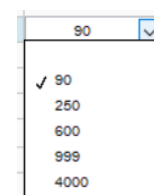
Fig. 10

### 6.1.6 Choosing the “Analysis Duration”

Using the **Analysis Duration**, it is possible to select different value for the duration of analysis. The selectable value depends on the sampling declared in “Sampling Configuration”.

### 6.1.7 Choosing the “Sampling Time”

Using the **Sampling Time**, it is possible to select different period of sampling. The selectable value depends on the sampling declared in “Sampling Configuration”.

## 6.2 Performing analysis

In order to start the analysis, press the START button consecutively 2 times (the former to turn on the display and the latter to start the analysis). Please refer to the RESPIROMETRIC Sensor operating instructions.

Before performing an analysis, all the steps described above must be carried out. Wait for the word “An” to appear on the sensor display, before starting another analysis.

The graph at the bottom of the display always refers to the analysis selected and not necessarily to the analysis in progress.

### 6.2.1 Stopping analysis

In order to stop the analysis, keep pressed the **SET+START** buttons at the same time. Please refer to the RESPIROMETRIC Sensor operating instructions.

Once the analysis has been stopped, the icon  will be showed on the left part of analysis. A stopped analysis can be validated by pressing the mouse right button and selecting the “Validate” option: the software will show the  icon and a message in the notes field.

### 6.2.2 Alarms

In order to see all the possible alarms, please refer to the RESPIROMETRIC Sensor operating instructions. RESPIROSoft™ gives on screen a short explanation of the problem.

### 6.3 Data evaluation

In order to evaluate the data obtained, the software offers the following possibilities:

- Configuring the data base
- Filtering the data base
- Exporting data

The indicator in the first column of the data base gives the status of each analysis in the series:

Indicator	Meaning
	Analysis ended correctly
	Analysis terminated manually with a reset by user.
	Analysis ended correctly but there are some notes (e.g. last value outside from range)

#### 6.3.1 Configuring the data base

The database consists of a set of data entered by the operator as well as data calculated by the software after analysis. Select **Modify** to modify the configuration of the database.

<b>Filter</b>	Filter data.
<b>Standard column</b>	View columns in standard configuration.
<b>Configure</b>	Change the order of columns (up and down arrows) or conceal columns (left and right arrows). Alternatively, press the right key of the mouse and select Configure (Fig. 11).
<b>Remove</b>	Delete one or more rows form the data base. Alternatively, press the right key of the mouse and select Remove.
<b>Select all</b>	Selecting all the rows in the database

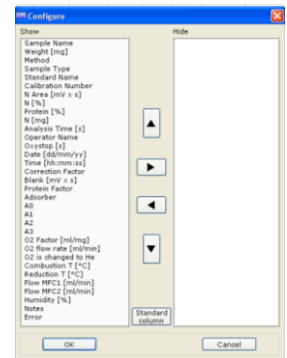
Fig. 11

#### 6.3.2 Filtering the database

The data shown can be filtered using the filtering function. Select **Modify/filters**.

Select filters as follows:

- 1- choose the fields to be filtered and the value required
- 2- press Save New in order to save the filter
- 3- Apply with Activate Filter or Deactivate Filter.

The filter status (“Filter ON” or “Filter OFF”) is shown at the bottom of the main window.

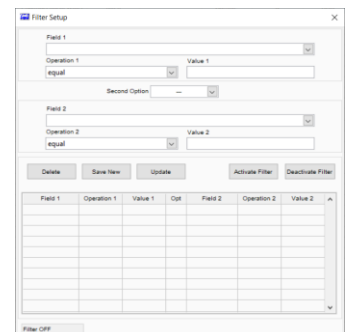
Fig. 12

## 7 Software installation requirements

<b>Operating system</b>	Windows 7 (SP 1), Windows 8 or Windows 10
<b>Hardware</b>	Core 2 Duo, RAM 4 GB
<b>Available space on the hard disk</b>	100 MB
<b>Recommended graphic resolution</b>	1024 x 768

## 8 Technical Data

Dimensions (WxHxD)	125x25x70 mm (4.9x1.0x2.7 in)
Weight	85 g (0.04 lb)
RESPIROMETRIC Sensor management	Up to 400 per Wireless DataBox™
Power supply	5 V – USB Cable
External power supply	100-240V±10% 50-60Hz 200mA
Power	0,225 W
Overvoltage category	II
Pollution degree	2
Max. altitude	2000 m
Max. humidity	85%
Admitted room temperature	Storing -25/+65° C Working +5/+40°C
Frequency WiFi	2.4GHz
Protocols WiFi	802.11 b/g/n

## 9 Spare Parts

10005112	USB cable A/Micro B 1.5m
40001862	Power supply 100-240V/5V no cord/plug
10003082	UE plug for power supply

## 10 Warranty

The unit is guaranteed against production defects for **25 months** from our invoice date. In accordance with this guarantee VELD Scientifica undertakes to repair any instruments resulting as faulty due to the quality of the materials used or poor workmanship. Faults arising due to inexpert handling/use or carelessness will not be replaced or repaired under. For more details please contact your Distributor.

### Exclusions:

The warranty will be considered null and void for faults resulting from:

- inexperience and carelessness of the customer.
- repairs, maintenance or replacements carried out by unauthorized third parties.
- use of non-original spare parts.

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## 11 Declaration of conformity

We, the manufacturer VELP Scientifica, under our responsibility declare that the product is manufactured in conformity with the following standards:

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EN 61010-1	Safety requirements for electrical equipment for measurement, control, and laboratory use Part 1: General requirements
EN 62311	Assessment of the compliance of low-power electronic and electrical equipment with the basic restrictions related to human exposure to electromagnetic fields (0 Hz to 300 GHz)
ETSI EN 301 489-1 V2.2.0	Electromagnetic compatibility (EMC) standard for radio equipment and services; Part 1: Common technical requirements
ETSI EN 301 489-17 V3.2.0	Electromagnetic Compatibility (EMC) standard for radio equipment and services; Part 17: Specific conditions for Broadband Data Transmission Systems
ETSI EN 300 328 V1.9.1	Wideband transmission systems; Data transmission equipment operating in the 2,4 GHz band

and satisfies the essential requirements of the following directives:

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2006/42/EC	Machinery directive
2014/53/EU	Radio Equipment directive (RED)
2015/863/EU (RoHS III)	Restriction of the use of certain hazardous substances
2012/19/EU (WEEE)	Waste of electric and electronic equipment

plus modifications / più modifiche / plus modifications / más sucesivas modificaciones / in der jeweils gültigen Fassung.

## 12 Declaration of conformity <sup>UK</sup> <sub>CA</sub>

We, the manufacturer VELP Scientifica, under our responsibility declare that the product is manufactured in conformity with the following regulations:

S.I. 2017/1206                      Radio Equipment Regulations 2017

according to the relevant designated standards:

EN 61010-1                      Safety requirements for electrical equipment for measurement, control, and laboratory use  
Part 1: General requirements

EN 62311                          Assessment of the compliance of low-power electronic and electrical equipment with the  
basic restrictions related to human exposure to electromagnetic fields (0 Hz to 300 GHz)

ETSI EN 301 489-1 V2.2.0      Electromagnetic compatibility (EMC) standard for radio equipment and services;  
Part 1: Common technical requirements

ETSI EN 301 489-17 V3.2.0    ElectroMagnetic Compatibility (EMC) standard for radio equipment and services;  
Part 17: Specific conditions for Broadband Data Transmission Systems

EN 61326-1                      Electrical equipment for measurement, control and laboratory use - EMC requirements -  
Part 1: General requirements

ETSI EN 300 328 V1.9.1        Wideband transmission systems; Data transmission equipment operating in the 2,4 GHz  
band

and satisfies the essential requirements of regulations:

S.I. 2008/1597                      Supply of Machinery (Safety) Regulations 2008

S.I. 2012/3032                      Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic  
Equipment Regulations 2012

S.I. 2013/3113                      Waste Electrical and Electronic Equipment Regulations 2013

plus modifications.

## Thank you for having chosen VELP!

Established in 1983, VELP is today one of the world's leading manufacturer of analytical instruments and laboratory equipment that has made an impact on the world-wide market with Italian products renowned for innovation, design and premium connectivity. VELP works according to **ISO 9001**, **ISO14001** and **ISO 45001** Quality System Certification.

Our instruments are manufactured in Italy according to the IEC 1010-1 and CE regulation.

Our product lines:

### Analytical instruments

Elemental Analyzers  
Digestion Units  
Distillation Units  
Solvent Extractors  
Fiber Analyzers  
Dietary Fiber Analyzers  
Oxidation Stability Reactor  
Consumables

### Laboratory Equipment

Magnetic Stirrers  
Heating Magnetic Stirrers  
Heating Plates  
Overhead stirrers  
Vortex mixers  
Homogenizers  
COD Thermoreactors  
BOD and Respirometers  
Cooled Incubators  
Flocculators  
Overhead Shakers  
Turbidimeter  
Radiation Detector  
Open Circulating Baths  
Pumps

[www.velp.com](http://www.velp.com)

VELP Scientifica Srl  
20865 Usmate (MB) ITALY  
Via Stazione, 16  
Tel. [+39 039 62 88 11](tel:+39039628811)  
Fax. [+39 039 62 88 120](tel:+390396288120)

*We respect the environment by printing our manuals on recycled paper.*

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